




**dentimarc**  
 Marco Millán González  
 Ca N°Andal, 53 08130 Santa Perpetua  
 Tel. 696 49 57 17 marco.dentimarc@gmail.com



**Revo-S™**  
**révolutionnaire !**



## Preparación y formación del canal

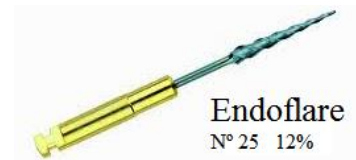

**SC1** 

(Shaper® & Cleaner 1)  
N°25 6%



**SC2** 

(Shaper® & Cleaner 2)  
N°25 4%



## Finalización del canal

**SU** 

(Shaper® Universal)  
N°25 6%




Protocole Protocol

**Revo-S® SC1 SC2 SU**

 Revo-S®

**1** SC1  2/3 LT-WL

**2** SC2  LT-WL


**3** SU  LT-WL

Vitesse 250 à 400 tr/min - Speed 250 to 400 rpm


Nous inventons, vous réussissez ! 


Protocole Protocol


**Revo-S® AS30 AS35 AS40**

 Revo-S®

Apex 30/100 → Apex 35/100 → Apex 40/100

**AS30**  LT-WL

**AS35**  LT-WL

**AS40**  LT-WL

Vitesse 250 à 400 tr/min - Speed 250 to 400 rpm